



Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model pledt-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control



\*Specifications subject to change without notice

Steridium Pty Ltd | Australia  
[www.steridium.com](http://www.steridium.com)  
[sales@steridium.com](mailto:sales@steridium.com)  
+61 7 55 931 352



Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model pledt-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control

Chamber capacity	2400 Litres
Standard dual chamber design	1 set of shelves per door
Chamber width	865 mm
Chamber depth	900 mm
Chamber height	1434 mm
Single chamber design	Optional
Chamber width	1809 mm
Chamber depth	900 mm
Chamber height	1434 mm
Chamber material	Stainless Steel 304L
External width	2214 mm
External depth	1067 mm
External height	2056 mm
External material	Durable epoxy powder-coated steel
Recommended clearance	75mm on sides and rear
Temperature range	4 to 50C
Temperature sensor	Pt100
Control accuracy	0.1 C
Temporal stability	0.2 C
Available shelf positions	30 (15 per door)
Shelves supplied	6 (3 per door)
Shelf dimensions	690 x 590 mm
Shelf spacing	87 mm
Shelf capacity	50 KG
Insulation type	Mineral Wool - asbestos free
Insulation thickness	50mm




Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model pledt-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control

Compressor size	1 HP
Refrigerant	CFC-Free R134a
Lighting source	Waterproof LED - lighting directly above entire shelf Mounted in lighting cassettes - height adjustable
Growth height	90 - 270mm height from shelf to light source
Growth area	4.6m <sup>2</sup> (6shelves)
Maximum intensity	150µmol/m <sup>2</sup> /sec
Lighting colours	Lighting cassettes comprised of any combination of: Warm White Cool White Red Green Blue RGB Combination—Variable spectrum

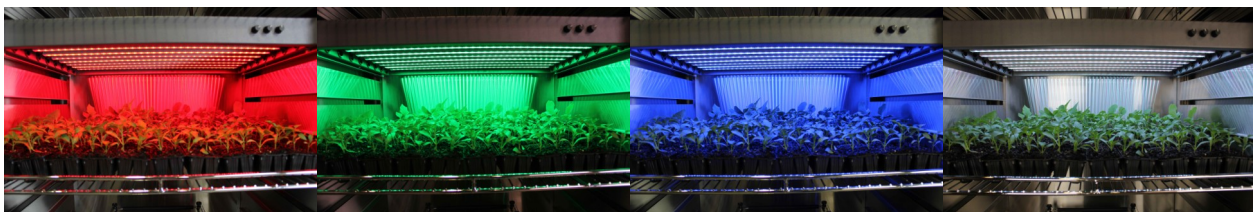



Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model plect-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control

Electrical supply	Single phase, 220 - 240V, 50 - 60 Hz (alternatives available on request)
Element rating	2 kW
Circulating fan	2
Controller type	Multifunction intelligent digital PID-microprocessor
Readout	Digital LCD readout of set and actual temperature
Temperature alarms	High   low   deviation
Programmable control	500 step programmable control of temperature, humidity and lighting
Country of origin	Australia




Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model pledt-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control

Humidity Range	20 - 90% RH
Humidity sensor	
Control accuracy	0.1 %
Temporal stability	0.5 %

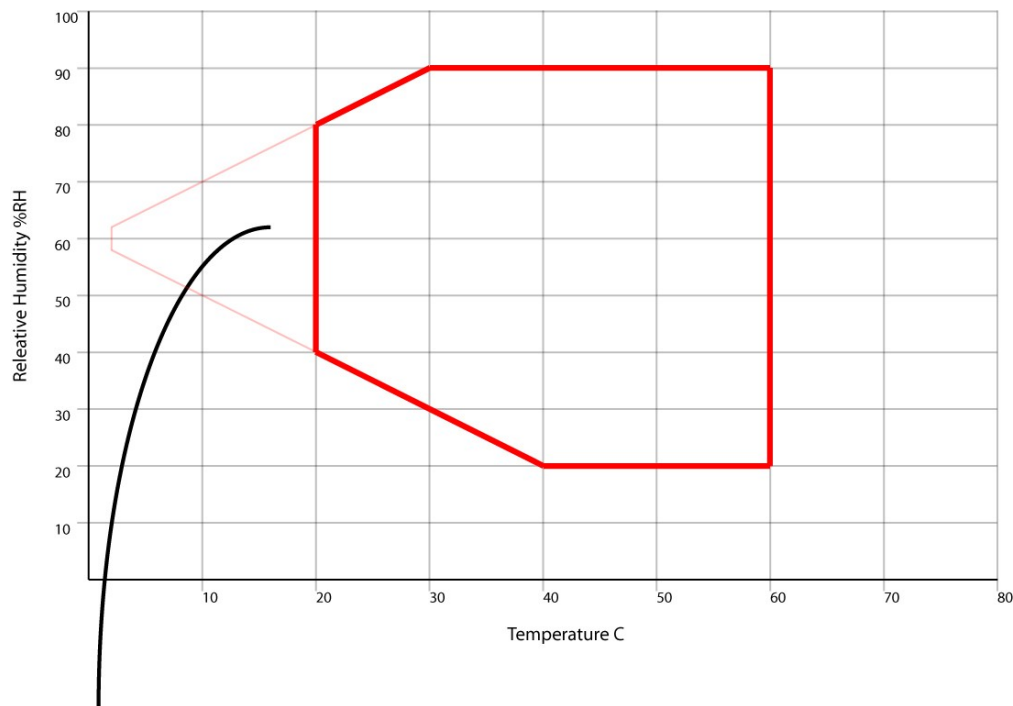

If the temperature setpoint is below 20 or above 80C, the RH control will be disabled and the chamber will adopt a natural value depending on the internal chamber conditions

This graph is a guide only  
Variations to the envelope can occur due to ambient conditions,  
the intensity of internal lighting, heat loads or humidity loads,  
and the nature of the product placed within the chamber.



Australia's premium manufacturer  
of medical and laboratory products

## Product Specification Sheet

Series Plant growth chambers  
LED Lighting  
Temperature and Humidity  
Model pledt-rh-2400

2400 litre LED plant growth chamber  
Temperature and humidity control

### Options / Accessories

- Single chamber design
- Additional shelves
- 8 Step programmable controller
- 500 Step programmable controller
- Inner viewing door - toughened glass
- Double glazed glass door
- Side access port - 50 mm
- Fan speed controller
- Internal power supply\*
- RS232 or RS485 comms module
- Ethernet comms module
- ISPM 15 compliant crating for export
- Calibration certificate
- Co2 Control (ambient to 3000ppm)

\*Maximum temperature limited to 50C